

ACCUderma  
BI-ENERGY



**Lifting in just  
one hour of treatment**



TERMO SALUD | MEDICAL  
& ESTHETIC

**ACCUderma bi-energy** is an effective solution that improves the **regeneration of collagen and elastin** fibres and achieves an **immediate lifting effect** by means of energetic and toning microcurrents.

It is a natural, non-invasive treatment that smooths the skin and improves the tone and luminosity. **ACCUderma bi-energy** acts directly on the connective tissues, thereby slowing down **the signs of aging**.

**Microcurrents** are a powerful and effective tool in the anti-aging battle as they reduce the appearance of fine lines and wrinkles, **improving the texture and appearance of the skin**.

## MODE OF ACTION


The **ACCUderma bi-energy** treatment applies low frequency, **low intensity variable microcurrents, emitted in short pulses**, in order to tighten the skin, reducing the signs of skin aging; such as wrinkles and flaccidity.


**ACCUderma bi-energy** activates the cellular system, due to an ionic sweep that causes soft, deep stimulation of the **facial tensor muscles**.

The equipment includes two types of microcurrents:


**Lift:** this step performs the function of preheating the cells by stimulating the superficial layers of the skin (dermis and epidermis) as a preliminary step prior to deeper stimulation.

**Tonify:** this step performs deep stimulation in the main tensor muscles.

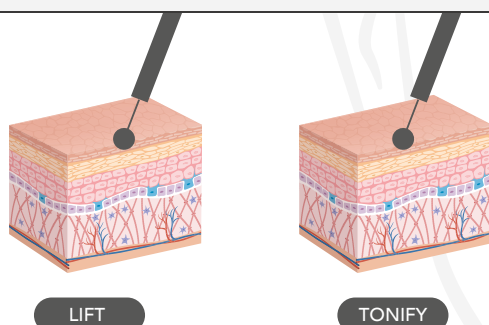

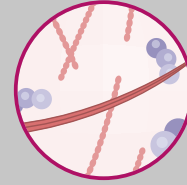


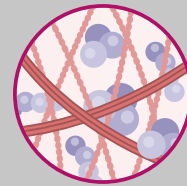

before  
**ACCUderma**  
BI-ENERGY



after  
**ACCUderma**  
BI-ENERGY

-  HYALURONIC ACID
-  COLLAGEN
-  ELASTIN

## BENEFITS

- Immediately oxygenates, repairs, regenerates, tones and hydrates the skin
- Increases the production of collagen and elastin
- Improves the texture and appearance of skin
- Improves skin microcirculation
- Activates the cellular system
- Stimulates the facial tensor muscles
- Drains accumulated toxins and fluids
- Relaxes the nervous system
- Visible results from the very first session
- Non-invasive treatment suitable for all skin types

## RESULTS

- Tissues are noticeably toned and defined
- Oxygenation of the skin increases and enhances luminosity
- The depth of wrinkles and expression lines decreases
- Skin gains firmness and muscle tone
- Decongests the skin

**ACCUderma**  
BI-ENERGY



**Lifting in just  
one hour of treatment**



Power supply:	Input tension: 100 -240 Vac, 50/60 Hz 0,6 A Output tension: 12 Vdc, 2,08A 25 W Max 6W
Maximum consumption:	Displays at 7 segments and rotary control wheel
Interface	IPX0
IP protection degree :	Electrical class I
Classification:	

## OUTPUT CHARACTERISTICS

### TONIFY

Output frequency: 1250Hz  
Intensity levels: 1 a 99  
Pulse duration: 400us  
Work cycle: 24ms  
Rest: 310ms  
Voltage: 60Vpp (en vacío)  
Treatment duration: 1 a 60min  
Electrical insulation class: Tipo BF

### LIFT


Output frequency: 185Hz  
Intensity levels: 1 a 99  
Pulse duration: 200us  
Work cycle: 500ms  
Rest: 1s  
Voltage: 60Vpp (empty)  
Treatment duration: 1 to 60 min  
Electrical insulation class: Type BF

Dimensions:	300 mm (h) x 260 mm (l) x 200 (w)
Weight:	3.9 Kg



Ref.: M97

TERMO SALUD | MEDICAL & ESTHETIC

 **Termosalud**  
Ataulfo Frieria Tarfe 8. 33211 Gijón. Asturias  
tel. +34 985 167 547  
email: [info@termosalud.com](mailto:info@termosalud.com)  
[termosalud.com](http://termosalud.com)

